BLACKROLL®

FASCIA GUN

Instruction Manual

EN | DE | FR | NL | IT | PL | CN

Safety Instructions

- + Do not use in sensitive areas or near organs.
- + If pain or any other contraindications occur, stop use immediately and consult a doctor.
- + People with reduced physical, sensory or mental capacities or anyone under the age of 18 must be supervised when using this device.
- + Do not put any electrically conductive objects into the slots of the device.
- + Do not insert the device into any of your body's natural orifices.
- + Never block the vibration of the device. Free vibration must always be possible.
- + Do not use as support or as a substitute for medical treatment without consulting your doctor.
- + With any use in the area of the head, neck and upper back, use caution and a targeted approach when applying vibration.
- + Do not immerse the device in water.
- + Keep liquids away from the ventilation openings and operating keys.
- + Keep the device away from sources of heat, children and pets.
- + Protect the device from moisture, dust, aggressive liquids and vapors.
- + If moisture does penetrate the device, switch it off immediately and if necessary disconnect the plug from the electrical outlet.
- + Do not remove any screws and do not attempt to dismantle the device.
- + Only pass the device on to third parties with the instructions for use.
- + Any warranty claims will be rendered invalid by cases of misuse.
- + Not to be used by diabetics.

Contraindications

Essentially there are hard indications that exclude the use of the device and contraindications in which a decision has to be made in each individual case whether the device should be used or not. Should you have one or more of the conditions listed below, please consult your doctor before using the device:

- + Redness
- + Swelling
- + Acute pain
- + Osteoporosis
- + Spinal disc damage
- + Taking anticoagulant (blood-thinning) medication
- + Open, inflamed or infected areas of skin
- + New fractures
- + Increased tendency to bleed, tendency to develop hematomas (bruises)
- + Tendency to develop embolisms
- + Cardiac pacemaker (PM)
- + Direct bone or eye contact
- + Fibromyalgia
- + Rheumatic disorders
- + Joint replacement
- + Tumor diseases
- + Pregnancy

Instructions for use

A

On/Off button & changing
massage levels

B

USB-C port

C


Status indicator: When switching on, the
battery status is shown. Each blinking LED at
C represents 25% battery level. During use, it shows
which of the four massage levels is currently selected.

D

Charging indicator: Lights up red during charging or when the
battery level is below 10%. Lights up green when the FASCIA GUN is
fully charged.

E

Interchangeable attachments

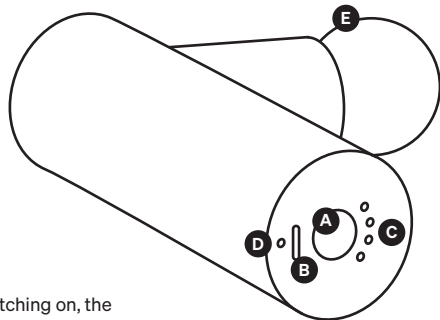1. Check all parts for damage.
2. Fully charge the battery using the USB-C port and the supplied cable before the first use. The device must not be used while charging. Charging time approx. 5h.
3. To switch the device on or off, press and hold the On/Off button on the bottom for 1.5 seconds. The LEDs on the status indicator will light up and display the battery level.
4. When switched on, you can change the massage levels by briefly pressing the On/Off button.

Massage levels & massage times

Adjust the level and duration of your massage as needed. You can adjust the massage mode and duration to suit your needs.

Level 1

Relax Mode

20 Hz; 1200r/min (3-4 minutes per body area)

Relax Mode uses a very low vibration frequency. It is ideal for relaxation and regeneration. Perfect for recovery or before going to bed. Always use low intensity pressure for this level.

Level 2

Post Exercise Mode

28.3 Hz; 1700r/min (2-3 minutes per body area)

The second level, Post Exercise Mode, is ideal for cooling down and preventing muscle soreness after an intense training session. The vibration can improve circulation and desensitise painful areas. This helps with your recovery.

Level 3

Pre Exercise Mode

45 Hz; 2700r/min (30 seconds - 2 minutes per body area)

The Pre Exercise Mode is designed to prepare you for exercise. This frequency improves circulation, tissue fluid flow and tissue quality. It means you gain increased mobility without decreasing your performance. The risk of injury can be reduced.

Level 4



Pro Mode

53.3 Hz; 3200r/min (30 - 60 seconds per body area)

Not for sensitive people. Pro Mode is very intense and can be used for large muscle groups such as the thighs or glutes. Use Pro Mode before training or competition.


Massage heads

The product comes with four different massage heads. You can easily change the attachments by pulling them off and putting a new one in.

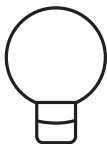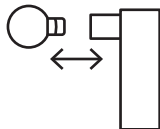Silicone Ball

Suitable for overall body treatments. Matches very well with our 1200r/min Relax Mode. The silicon surface ensures enough grip for a targeted massage.

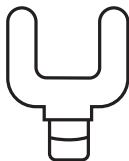Silicone Fork

Perfect for the massage of parallel muscles, left and right of the Achilles tendon or the upper area of the trapezius muscle.

Silicone Trigger

The right attachment for selective trigger point massage with vibration. Helps to release stiff tissues or fascia.

Metal Flat

Can be used for overall body treatment. Gliding metal material enables stroking movements on textile or skin. Also suitable to use with massage oil.

Cleaning

- + Do not clean when product is charging.
- + Do not use water to clean the product. Instead use a wet cloth.
- + Do not use any strong detergents or solvents.
- + Attachments and device can be cleaned with disinfectant.

Technical specifications

Battery capacity in mAh:	2550
Maximum charging current in A:	2
Maximum charging voltage in V:	5
Type of battery:	Li-Ion
Can the battery be replaced?	no
Type of charging connection:	USB-C
Is a charging cable supplied?	yes
Power adapter included?	no, please use 5V-2A power adapter
Protection class:	IP22
Standby:	15h in RELAX MODE (without load)
Material:	Aviation aluminum alloy
Stalling force in KGF:	22
Net weight in g:	550
Dimensions in mm:	150*92.5*48
Amplitude in mm:	6

Further instructions for use

Self-massage with the FASCIA GUN is very intuitive. This means that you yourself know best where your tissues need attention.

It is generally best to glide the FASCIA GUN gently over the tissue. Stay on sensitive hotspots/trigger points for 30 seconds or until tensions ease.

You can use the FASCIA GUN to passively massage tissue by simply applying pressure to your tissue. Or you can add active mobilisation by moving adjacent body parts at the same time.

For example, while massaging your calf, you can simultaneously bend or stretch your ankle for an additional mobilisation effect. This combined movement creates shear forces in the tissues and is more effective in loosening tensions.

Avoid bony areas. Self-massage is designed for the active musculoskeletal system (muscles, tendons and fascia) and not for the passive musculoskeletal system (bones, joints, intervertebral discs). We recommend a low vibration level if you are working close to a bone.

Areas that are difficult to reach (e.g. back and neck) are best massaged by a partner.

You can vary the intensity of the massage by adjusting the angle at which you work your muscles. When applied vertically, the gentle strokes reach deeper into the tissue. Using the FASCIA GUN at an angle reduces the intensity and is especially suitable for sensitive parts of the body.

You can use the FASCIA GUN daily for self-massage. We recommend giving yourself a day off if your tissues feel a little stressed the next day. The best results are achieved with regular use.

Note: The information given is only a recommendation. Use the product in a way that is comfortable for you. If you are unsure about how to use the product correctly, please contact your doctor or physiotherapist.

Storage

Please store the gun in a dry place away from heat, children and pets. Use the packaging or the FASCIA GUN CASE while you are not using the gun.

Disposal

This product is subject to the European Directive 2012/19/EU (WEEE) and must not be disposed of with household waste. Electrical and electronic equipment must be disposed of at designated collection points or take-back systems. Information on collection points in your area can be obtained from your local municipal waste disposal authority.

The lithium-ion battery contained in this device is subject to the European Directive 2006/66/EC and the Battery Act. Lithium-ion batteries must not be disposed of in household waste. You are obliged to return the battery to suitable collection points or to the retailer free of charge. You can find information about return points from your retailer or from municipal disposal facilities.

Important note: Lithium-ion batteries can pose a risk to the environment and health if disposed of improperly. Batteries must not be damaged, heated or thrown into a fire.

Disclaimer

The unit is used for exercise at your own risk. It is recommended that you seek advice from a physiotherapist, manual therapist, masseur or medical specialist if you experience severe pain before or after self-massage or recurrent injuries.

Warranty

The seller grants the buyer, if he is a consumer, the statutory warranty in accordance with the legal requirements for a period of 2 years from the purchase of the item. The Seller's warranty obligations shall initially be limited to supplementary performance, whereby the Seller shall be free to choose whether to deliver a defect-free item or to repair the defective item. If subsequent performance is unsuccessful, the buyer shall be entitled to the further statutory claims. The seller rejects warranty claims for damage caused by accident, obviously incorrect use, misuse or improper transport; which occurred as a result of repairs carried out by unauthorised personnel; which occurred as a result of use and handling not in accordance with the operating instructions; which occurred as a result of improper handling of the purchased item. Furthermore, claims for damages by the purchaser are excluded. Excluded from this are claims for damages by the buyer arising from injury to life, limb or health as well as liability for other damages based on an intentional or grossly negligent breach of duty by the seller, his legal representatives or vicarious agents. Also unaffected is the liability for breaches of duty, the fulfilment of which makes the proper execution of the contract possible in the first place and on the observance of which the buyer may regularly rely. In the event of a slightly negligent breach of such contractual obligations, the Seller shall only be liable for the foreseeable damage typical for the contract, unless it is a matter of claims for damages by the Customer arising from injury to life, limb or health.

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BLACKROLL AG
WSVK Oederan GmbH